

# STAIR TRAC

## PORTABLE WHEELCHAIR LIFT



The Garaventa Stair-Trac is an easy to use, portable wheelchair lift designed to attach to most manual wheelchairs. The powerful motor allows a small attendant to safely transport a larger passenger up or down stairs. The Stair-trac is suitable for indoor or outdoor use, in public buildings or private homes. It is an affordable and reliable accessibility solution.



# STAIR-TRAC

## FEATURES

- Attaches to most standard manual wheelchairs with push handles
- Robust battery - fully charged it can carry 150 kg (330 lbs) up 30 flights of stairs
- Works on virtually any indoor or outdoor stairway
- Stable and self supporting
- Durable steel reinforced rubber tracks grip stairs securely
- Detachable support handle for compact storage and easy transport
- Quick release safety belt
- Key switch controls unauthorized use
- Rechargeable batteries and charger

## BENEFITS

- An economical alternative to stairway lifts and elevators
- May be used on several different stairways
- Provides accessibility where power is unavailable
- Stores easily when not in use
- Requires minimal maintenance
- Available immediately
- 2 year standard warranty
- Available extended warranty\* (3 additional years)  
\*USA/Canada only



The optional vehicle loading ramp allows the Stair-Trac to be easily and safely loaded and unloaded from vehicles using the drive motor.

## SPECIFICATIONS

### Dimensions:

Height: 903 mm (35 5/8")  
 Length: 1453 mm (57 1/4")  
 Width: 644 mm (25 3/8")

**Weight:** 59 kg (130 lbs)

**Load Capacity:** 150 kg (330 lbs)

### Speed:

Up: 6.1 meters (20 ft) per minute  
 Down: 10.7 meters (35 ft) per minute

### Operating Controls:

Constant pressure drive control switch

**Motor:** 12V DC

**Batteries:** 2 x 6V - 20A hours


**Maximum Charge Rate:** 4 amp

**Optimum Stair Angle:** 35°

**Typical Battery Charge Capacity:**  
650 stairs

Questions?

Call our Product Information Hotline at:  
1.800.663.6556

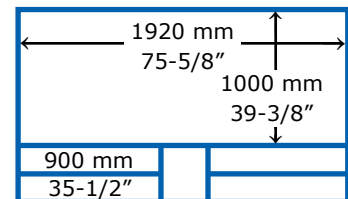

### Required Clearances:

This diagram (above) shows typical clearances required on turnback stairways. Dimensions vary slightly with stairway width and operator experience.



Authorized Garaventa Lift Representative

### Contact us today

Phone: +1 604 594 0422  
 Toll Free: 1 800 663 6556  
 Web site: [www.garaventalift.com](http://www.garaventalift.com)

